DPE-02

Dynamic Pressure Exciter







Advantages

- Realization of all automatic calibrations according to own test regime (up to 1 shock/s)
- User-friendly and efficient
- Upgradeable to a shock acceleration calibration device

Application

 Generates pressure shock impulses with low pulse duration (<2 ms)

Range of use

- Dynamic calibration of pressure sensors
- Dynamic tests of pressure sensors

Features

- Pressure amplitudes from **2 MPa to 420 Mpa** with pulse duration less than 2 ms
- Create test descriptions in which a calibration run is performed fully automated using the SPEKTRA CS18 Software

DPE-02

Dynamic Pressure Exciter




| Features | SPEKTRA Shock pressure device |
|-----------------------|------------------------------------|
| amplitude | 2 420 MPa |
| pulse duration | <2 ms (depends on projectile mass) |
| operation | pneumatic, fully automated |
| Impulse form | half sine |
| pressure transmission | hydraulic (Oil) |
| upgradeable | shock acceleration |

Schematic

Pressure chamber DUT Piston Projectile Barrel Linear Stage Pressurised air

Pulse form

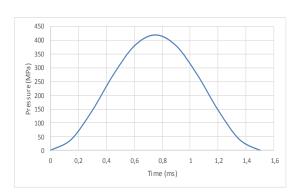

Subject to change without notice